ConeGuard™ Thread Protector


Extend the life of your ICP-MS Cones

When cleaning cones that have a screw thread, it is important that the thread is not contacted by any corrosive solution. If the thread gets corroded, the cone may not seal correctly or it may bond to the base and be difficult to remove. And with Pt cones, the thread is likely to wear out before the Pt insert.

It is also improtant to keep the thread in good condition to prevent cross-threading and potential damage to the cone housing.

The Glass Expansion ConeGuard Thread Protector seals the thread and protects it from corrosion during the cleaning process. The ConeGuard is available for a range of ICP-MS Models.



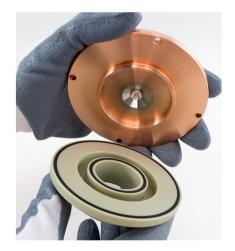



For more information visit: www.geicp.com or contact us at enquiries@geicp.com



6 Central Boulevard Port Melbourne VIC 3207 Australia

+61 3 9320 1111 enquiries@geicp.com

Americas

31 Jonathan Bourne Drive Unit 7, Pocasset, MA 02559 USA

508 563 1800 geusa@geicp.com

Europe

Friedenbachstrasse 9 35781 Weilburg Germany +49 6471 3778517

+49 6471 3778517 gegmbh@geicp.com